

hyperClone® Human Fresh Blood
NK Activation/Expansion Kit

Kit03



人新鲜血样来源NK高效版激活扩增试剂盒

NK > 产品性能




CD3⁻CD56⁺ NK细胞纯度 > 98% , CD3⁺ < 1% , 满足自体或异体使用

50ML外周血或脐带血, 扩增千亿级别NK细胞数量, 大量现货型NK产品


可持续培养超过28天, 细胞活率高于90%

满足NK细胞新药申报需求

NK > 数据展示

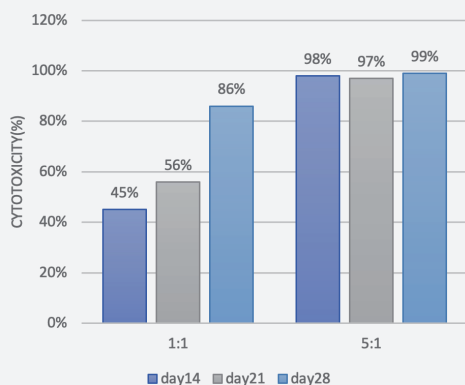

细胞杀伤结果

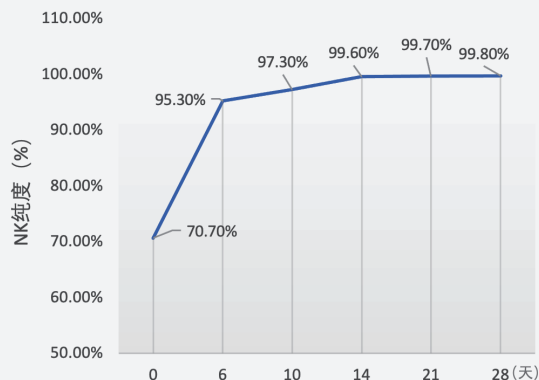

细胞纯度

Time (days)	Cell Survival Rate (%)
0	99.50%
6	93.40%
10	93.10%
14	96.90%
21	96.80%
28	96.30%

NK细胞扩增数量和扩增倍数

Time Point	扩增细胞数 (Expanded Cell Count)	扩增倍数 (Expansion Fold)
DAY0	8.98E+06	1
DAY6	2.31E+08	26
DAY8	6.90E+08	77
DAY9	2.60E+09	296
DAY10	4.47E+09	496
DAY12	8.64E+09	962
DAY14	1.84E+10	2043
DAY16	3.38E+10	3743
DAY18	6.67E+10	7423
DAY21	1.52E+11	1692

SSC-H :: SSC-Height

Lymphocytes
95.2

FSC-H :: FSC-Height

苏州华辰生物科技有限公司

Suzhou Huachen Biotechnology Co., Ltd

电话：400-965-9800

网址: www.huachenbiotech.com

地址：苏州工业园区星湖街328号创意产业园A3-504-1单元

